**Microbiological Incubator
VWR INCU-Line IL 10
VWR INCU-Line IL 23
VWR INCU-Line IL 56**

INSTRUCTION MANUAL

European Catalogue Numbers:

**390-0384 IL 10 digital, transparent door window
390-0385 IL 10 digital, opaque door window
390-0482 IL 23 digital, transparent door window
390-0999 IL 56 digital, transparent door window**

Version: 11

Issued: 22.11.2018

Legal Address of Manufacturer

Europe

VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
B-3001 Leuven
+ 32 16 385011
<http://be.vwr.com>

Country of origin: SLOVENIA

Table of Contents

Warning	4
Safety information	4
Package Contents	4
Unpacking.....	5
Installation.....	5
Intended Use	5
Symbols and conventions	5
Product specifications	6
Overview.....	6
Description of front panel	6
Instructions for use	7
Getting started	7
Operation	7
Temperature re-adjustment.....	7
Troubleshooting	8
Cleaning and maintenance	8
Accessories	8
Technical Service	9
Warranty	9
Compliance with local laws and regulations.....	9
Equipment disposal	9
Appendix.....	10
Your Distributor.....	12

Warning

Configuring and operating this equipment requires **precise knowledge** and **observance** of the instructions in this **user manual**. The equipment is designed only for its intended application.

Important safety symbols used in this manual:

Warning!

This is a warning of risk situations and dangers. Failure to observe this warning could be life-threatening. These warnings must be heeded.

This equipment is designed for 230V (50/60Hz).

Make sure it is connected to the correct power supply.

Disconnect the power cord before carrying out maintenance or cleaning work.

Safety Information

The responsibility for operating the equipment lies exclusively with the owner or user if said equipment is incorrectly serviced, maintained or altered by persons not employed by an authorized dealer, or if the equipment is used in a manner contrary to its specified purpose.

The equipment must be maintained and operated in accordance with this user manual.

Ensure that this user manual is accessible at all times and has been read and understood.

To avoid electric shock do not spill any liquids on the equipment.

Package Contents

IL 10:

Description	Quantity
Mini Incubator IL 10	1
Stainless steel shelf	1
Power cord	1

IL 23:

Description	Quantity
Incubator IL 23	1
Stainless steel shelf	2
Power cord	1

IL 56:

Description	Quantity
Incubator IL 56	1
Stainless steel shelf	2
Power cord	1

Unpacking

Carefully inspect the unit and accessories thoroughly upon receipt for any damage that may have occurred during shipping. Severe damage to the packaging container may indicate damage to the unit. If you suspect damage to the unit, immediately file a claim with the carrier in accordance with their instructions before contacting VWR. Also verify that all accessories are included and that the unit is in good working order before discarding the shipping packaging.

Installation

Your satisfaction and safety require a complete understanding of this unit. Read the instructions thoroughly and be sure all operators are given adequate training before attempting to put the unit in use.

NOTE: This equipment must be used only for its intended application; any alterations or modifications will void your warranty and may cause injury.

Intended use

VWR INCU-Line IL 10 & IL 23 & IL 56 are compact general purpose convection incubators capable of uniformly heating samples up to a maximum of 70 °C. The units have a door window area so samples can be observed without opening the door, thus preventing heat loss.

For research use only. Not intended for any animal or human therapeutic or diagnostic use.

Symbols and conventions

The following chart is an illustrated glossary of the symbols that are used in this manual.


	CAUTION This symbol indicates a potential risk and alerts you to proceed with caution.
	CAUTION This symbol indicates the presence of high voltage and warns the user to proceed with caution.
	CAUTION This symbol indicates risks associated with hot surfaces.

Product Specifications


Specifications	IL 10	IL 23	IL 56
Electrical	230V, 50/60 Hz, 0.3 A	230V, 50/60 Hz, 0.7 A	230V, 50/60 Hz, 0.96 A
Power consumption	70 W	160 W	220 W
Fuses	2 x 0.5 AT 250 V	2 x 1 AT 250 V	2 x 2 AT 250 V
Controls	PID control	PID control	PID control
Temperature range	Ambient +5 °C to +70 °C	Ambient +5 °C to +70 °C	Ambient +5 °C to +70 °C
Temperature resolution (display)	0.1 °C	0.1 °C	0.1 °C
Temperature accuracy at 37 °C	± 0.5 °C	± 0.5 °C	± 0.5 °C
Temperature fluctuation (time) at 37 °C	± 0.5 °C	± 0.5 °C	± 0.5 °C
Temperature variation (spatial) at 37 °C	± 0.5 °C	± 0.5 °C	± 1 °C
Environment temperature	+5 °C to +40 °C	+5 °C to +40 °C	+5 °C to +40 °C
Interior volume	10 l	23 l	56 l
Maximum load per shelf	3 kg	3 kg	10 kg
Maximum load per chamber	6 kg	9 kg	30 kg
Exterior dimensions W x D x H	285 x 280 x 335 mm	340 x 360 x 400 mm	453 x 420 x 546 mm
Interior dimensions W x D x H	230 x 220 x 200 mm	290 x 290 x 270 mm	390 x 360 x 400 mm
Weight	8.3 kg	12.8 kg	23.6 kg

Overview

Access port (Ø 7 mm) for external temperature probe

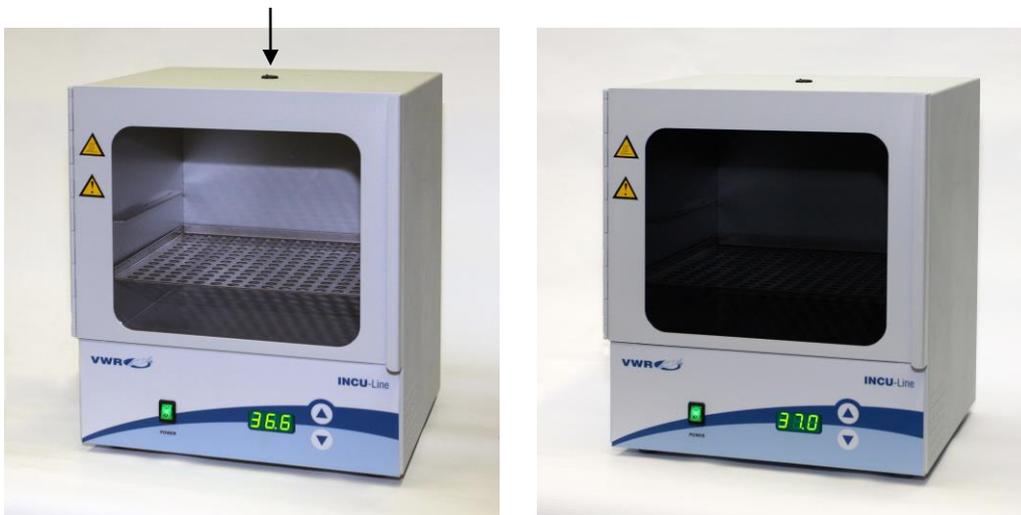Description of front panel

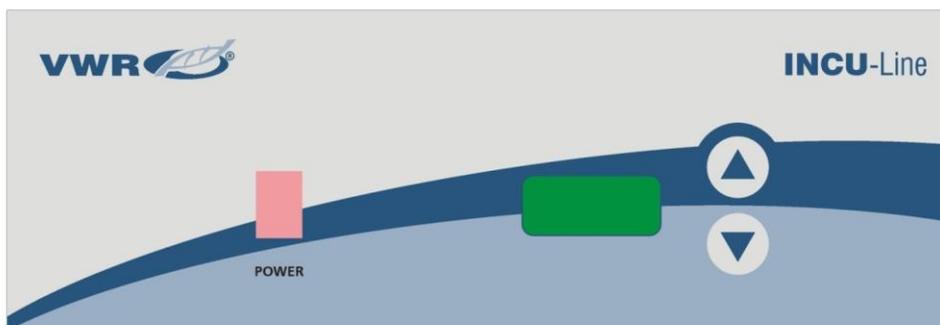POWER - power switch: For turning the incubator ON and OFF. It lights when incubator is switched ON.

▲ - UP key: For increasing the temperature.

▼ - DOWN key: For decreasing the temperature.

Display: For displaying the temperature (dot blinks when you set the temperature).

Instructions for use

Getting Started

Check the data plate for voltage, cycle, wattage and ampere requirements. If matched to your power source, plug the power cord into a grounded outlet. Voltage should not vary more than 10% from the data plate rating.

In selecting a location, consider all conditions which might affect performance, such as heat from radiators, autoclaves, etc. Avoid direct sun, fast-moving air currents, heating/cooling ducts and high traffic areas. Allow a minimum of 5 cm between the unit and walls or partitions which might obstruct free airflow.

Operation

Do not use this incubator in explosive or flammable environments. Do not heat or incubate flammable, explosive or highly reactive materials in this incubator as serious injury may result.

Place the unit on a level surface and plug the unit into a properly grounded outlet of appropriate voltage.

Install the shelf at the desired height, with the shelf edges up. If you intend to place the samples on the bottom of the chamber, install the shelf so that it sits with the edges down on the chamber bottom, thus creating an air gap between the chamber bottom and your samples.

Place the samples in the chamber and close the door. Make sure the door is securely closed.

Turn the power switch to the ON position. The switch should illuminate. On display all "888" are displayed at first, then automatically detected power supply frequency "F50" or "F60", then software version "VX.X", and finally the current temperature in the chamber. Incubator starts to heat.

Set desired temperature with UP or DOWN key. You must hold UP or DOWN key more than 2 sec. DOT on DISPLAY blink. When DOT start blink (dot blink 3 sec) set temperature on desired value between room temperature +5 °C to +70 °C. When desired temperature is set, display will after 3 sec automatically show momentary temperature.

You can always look at SET temperature when hold UP or DOWN key more than 2 sec.

Temperature re-adjustment

Incubator has software for temperature re-adjustment. You must first measure temperature in middle of chamber or in sample for minimum 2 hours. You must measure temperature with a certificated reference temperature measuring instrument with precision 0.1 °C or more. After 2 hours or more read temperature on thermometer and compare it with temperature on LED display. Difference from thermometer and display is value, which you will enter in software.

Example 1: Reference temperature is 37.9 °C, temperature on display is 37 °C. Difference is $37.9 - 37 = 0.9$
Value is 0.9 (value which you enter in software).

Example 2: Reference temperature is 36.2 °C, temperature on display is 37 °C. Difference is $36.2 - 37 = -0.8$
Value is -0.8 (value which you enter in software).

Procedure for temperature re-adjustment:

Press UP and DOWN key at the same time for more than 5 seconds - on display appears "COr" (correction) and DOT blink.

With UP or DOWN key define value (when dot blink). Pay attention on sign of value (plus or minus). After 3 sec when you will define temperature display will automatically show temporary temperature.

Troubleshooting

Review the information in the table below to troubleshoot operating problems.

Problem on display	Cause	Solution
Er1	Temperature probe fault, electronic board fault or incorrect setting of parameters.	Contact your distributor.
Er2	Momentary temperature in the incubator is at least 5 °C higher than SET temperature, because the incubator has been working at higher SET temperature recently.	Check if SET temperature is for 5 °C, or more, lower than momentary temperature in the incubator. If so, allow incubator to cool down. Turn the power switch off and on, then try again.
	Momentary temperature in the incubator is at least 5 °C higher than SET temperature, due to temperature probe fault, electronic board fault or incorrect setting of parameters.	Contact your distributor.
Er3	Momentary temperature does not reach SET temperature in 2 hours, because the door is opened, you open the door too frequently or the environment temperature is too low.	Check that the door is securely closed, don't open the door too frequently during operation, check that environment temperature is above minimum. Turn the power switch off and on, then try again.
	Momentary temperature does not reach SET temperature in 2 hours, due to thermal fuse cut off, heater fault, electronic board fault or incorrect setting of parameters.	Contact your distributor.

Cleaning and maintenance

Prior to any cleaning or maintenance of the unit, disconnect the power cord from the wall outlet.

No routine maintenance is required on the unit.

Clean the unit internally and externally with a damp cloth and mild soap. Do not immerse the unit or spill liquids into or on the unit. Use of alcohol or ammonia-based cleaners on the window may cause the window to crack. Allow the incubator to dry completely inside and out before reconnecting the power and using.

Accessories

IL 10:

Description	Quantity	Cat. No.
Extra shelf, 226 x 215 mm, stainless steel	1	390-0387
Tube rack, 2 x 6 positions, Ø 32 mm	1	390-0386

IL 23:

Description	Quantity	Cat. No.
Extra shelf, 286 x 285 mm, stainless steel	1	390-0483
Tube rack, 2 x 6 positions, Ø 32 mm	1	390-0386

IL 56:

Description	Quantity	Cat. No.
Extra shelf, 385 x 345 mm, stainless steel	1	390-1016
Tube rack, 2 x 6 positions, Ø 32 mm	1	390-0386

Technical service

Web Resources

Visit the VWR's website at www.vwr.com for:

- Complete technical service contact information.
- Access to VWR's Online Catalogue, and information about accessories and related products.
- Additional product information and special offers.

Contact us

For information or technical assistance contact your local VWR representative or visit www.vwr.com.

Warranty

VWR International warrants that this product will be free from defects in material and workmanship for a period of two (2) years from date of delivery. If a defect is present, VWR will, at its option and cost, repair, replace, or refund the purchase price of this product to the customer, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear. If the required maintenance and inspection services are not performed according to the manuals and any local regulations, such warranty turns invalid, except to the extent, the defect of the product is not due to such non-performance.

Items being returned must be insured by the customer against possible damage or loss. This warranty shall be limited to the aforementioned remedies. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

Compliance with local laws and regulations

The customer is responsible for applying for and obtaining the necessary regulatory approvals or other authorizations necessary to run or use the Product in its local environment. VWR will not be held liable for any related omission or for not obtaining the required approval or authorization, unless any refusal is due to a defect of the product.

Equipment disposal

This equipment is marked with the crossed out wheeled bin symbol to indicate that this equipment must not be disposed of with unsorted waste.

Instead it's your responsibility to correctly dispose of your equipment at lifecycle-end by handling it over to an authorized facility for separate collection and recycling. It's also your responsibility to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect from health hazards the persons involved in the disposal and recycling of the equipment.

For more information about where you can drop off your waste of equipment, please contact your local dealer from whom you originally purchased this equipment.

By doing so, you will help to conserve natural and environmental resources and you will ensure that your equipment is recycled in a manner that protects human health.

Thank you

VWR International bvba
Haasrode Researchpark 2020, Geldenaaksebaan 464, B-3001 Leuven
<http://be.vwr.com>

EU Declaration of Conformity

Manufacturer: VWR International bvba, Haasrode Researchpark 2020, Geldenaaksebaan 464, B-3001 Leuven

The undersigned states that this declaration is issued under the sole responsibility of the manufacturer, VWR International bvba, and that the object(s) of the declaration described in annex 1 are in conformity with the relevant Union harmonization legislation.

In the event of unauthorized modification of any of the products listed in the annex 1, this declaration becomes invalid.

Object Name: VWR Incubators INCU-Line® IL

Relevant EC Directives:	Low Voltage Directive	2014/35/EU
	EMC Directive	2014/30/EU
	RoHS	2011/65/EU

Harmonized Standards:	EN 55011:2009
	EN 55014-1:2006
	EN 55014-2:2015
	EN 61000-3-2:2014
	EN 61000-3-3:2013
	EN 61000-6-1:2007
	EN 61326-1:2013
	EN 61010-1:2010
	EN 61010-2-010:2014
EN 50581:2012	

Document Number: 28429.2017.01

Place & date of issue: Leuven, 10/04/2017

Signed for and on behalf of VWR International bvba:

Yves Van Damme
European Regulatory Affairs Manager

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India

VWR Lab Products Private Limited
No.139, BDA Industrial Suburb,
6th Main, Tumkur Road, Peenya Post,
Bangalore, India – 560058
Tel.: +91-80-28078400
Email: vwr_india@vwr.com

Ireland / Northern Ireland

VWR International Ltd / VWR
International (Northern Ireland) Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: +353 01 88 22 222
Email sales.ie@vwr.com

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VWR International GmbH
VWR International China Co., Ltd
VWR Singapore Pte Ltd